By Bob Page, AAS, NREMT-P, CCEMT-P, I/C

Table of Contents	Page(s)
Introduction and Objectives	2-3
Module One: Review of Anatomy and Physiolo	рду
Airway Review	4-7
The Physiology of Respiration and Perfusion	8-14
Module Two: The Technology of Capnography	,
How EtCO ₂ is measured	15
Review of Technology	16
Module Three: Clinical Applications of Capnogra	aphy
Waveform Interpretation	17-26
ETT Conformation	27
Closed Head Injuries	27
Ventilatory assessment	27
Perfusion trending	27
Module Four: Practice Cases and Evaluations	28-34
Closing Remarks	35

INTRODUCTION

Capnography is a noninvasive method for monitoring the level of carbon dioxide in exhaled breath ($EtCO_2$), to assess a patient's ventilatory status. A true capnogram produces an $EtCO_2$ value as well as a waveform, or capnogram. On Critical Care transports, capnograms are useful for monitoring ventilator status, warning of airway leaks and ventilator circuit disconnections. Capnography is also useful for ensuring proper endotracheal tube placement. Capnography also helps clinicians diagnose specific medical conditions, make treatment decisions, and assess efficacy of code efforts and predict outcome.

Capnography offers numerous clinical uses, but technical limitations have prevented EMS personnel from embracing its use outside the operating room. Today, technological advances have made it possible for these devices to be used in the demanding setting of the prehospital environment.

OBJECTIVES for the Session:

By the end of this session, you will be able to:

Describe the structure and function of the upper and lower airways.

Describe the mechanics and science of ventilation and respiration.

Describe the basic physiology of perfusion.

Describe the relationship between ventilation and perfusion.

Describe the principles behind CO_2 measurement.

Describe the various methods of EtCO₂ measurement including quantitative and qualitative capnometry and capnography.

Describe the technology of EtCO₂ measurement including mainstream, sidestream and microstream sampling.


Identify the components of a normal capnogram waveform.

Identify abnormal capnogram waveforms as related to various airway, breathing and circulation problems.

Discuss the various clinical applications of capnography in the field.

Given various cases, discuss the role of capnography in identifying the problem and in the management of the patient.

MODULE ONE: REVIEW OF AIRWAY ANATOMY AND PHYSIOLOGY

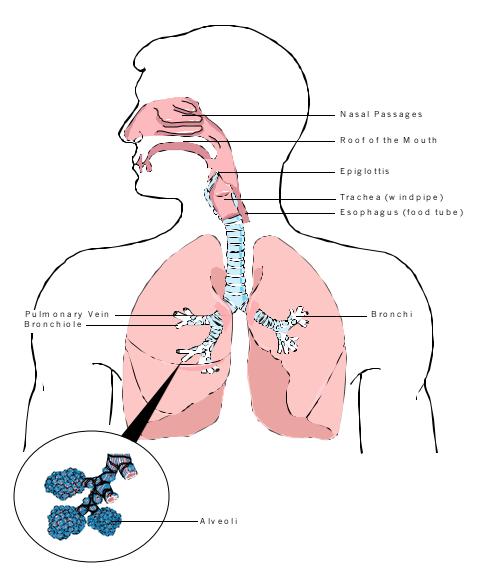

Two divisions to the airway:


- 1. Upper Airway
- 2. Lower Airway

Airway Anatomy Review

Upper Airway	
Nasopharynx	
Oropharnyx	- Nasopharynx - Oropharynx
Palates Hard and soft	Laryngopharymx
Larynx Epiglottis Vocal cords	


Upper Airway Physiology: PEEP or Positive End Expiratory Pressure can be defined as the pressure against which exhalation occurs. The purpose of PEEP is to prevent alveolar collapse. The structure and path way of the upper airways provide for a natural or "physiological PEEP."

The primary roles of the upper airway are to		
, and	the air entering	
the lungs.		

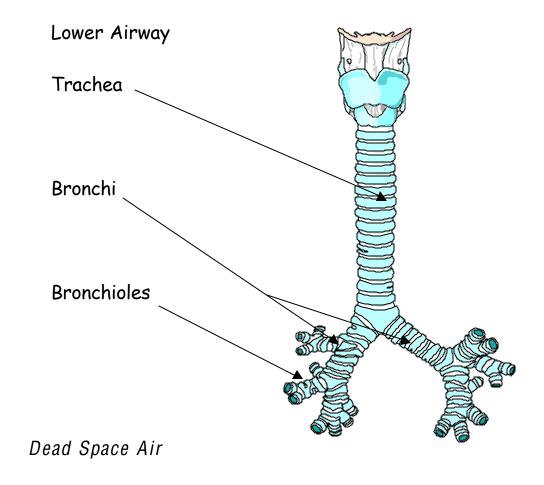


The lower airway is comprised of the trachea, bronchi and then 25 divisions of the bronchial tree terminating at the respiratory bronchioles and the alveoli.

The areas from the bronchioles to the nose comprise the total dead space air. This is air that is not exchanged.

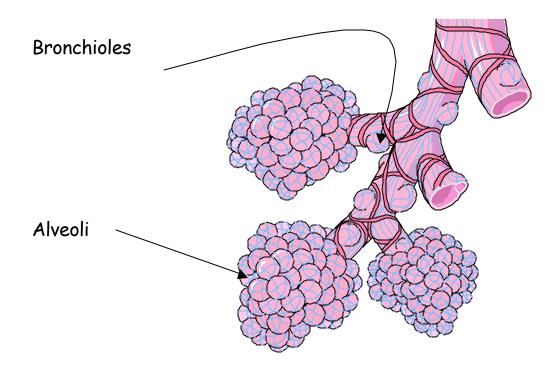

The alveoli are tiny air sacs where gas exchange occurs. O_2 and CO_2 are exchanges at the capillary-alveolar membrane.

THE PHYSIOLOGY OF VENTILATION

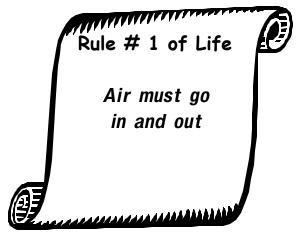

Function of ventilation:

Ventilation is the movement of air.

Designed to eliminate CO₂ and take in O₂

How air moves (and it MUST move)

Chemoreceptors in the medulla sense elevated levels of CO_2 or lowered pH, triggering ventilation. Known as hypercarbic drive.

Diaphragm contracts and moves downward.

Intercostal muscles spread chest wall out increasing the volume inside the chest.

Differences in pressure inside the chest and outside causes air to move into the lungs.

Hypoxic drive (low O_2 levels) is secondary drive.

VOLUME CAPACITIES

Tidal Volume (Vt): The amount of air moved in one breath Typically 500 cc in an adult at rest

Anatomical Dead Space (Vd): Air not available for gas exchange

About 150 cc

Alveolar volume (Va); Air that is available for gas exchange About 350 cc (Vt - Vd = Va)

Anything that affects the tidal volume only affects the alveolar volume.

Factors affecting the tidal volume

Hyperventilation

Fast breathing (tachypnea) doesn't necessarily increase tidal volume

Anxiety, head injuries, diabetic emergencies, PE, AMI, and others

Hypoventilation


Slow breathing (bradypnea) does not necessarily decrease tidal volume.

Causes include CNS disorders, narcotic use and others.

Street Wisdom: An increase or decrease in tidal volume is at the expense or benefit of alveolar air.


Respiration is the exchange of gases

Alveolar respiration occurs between the _____ and ____ in the lungs.

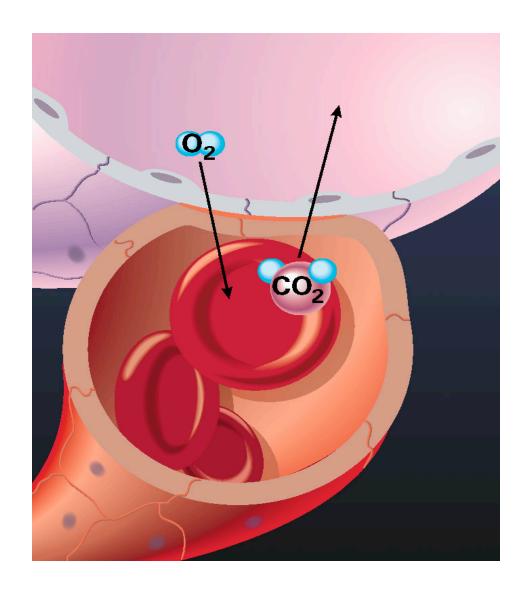


Cellular respiration occurs between the _____ in the body and the _____.

When there is a difference in partial pressure between the two containers, gas will move from the area of greater concentrations to the area of lower concentration. a.k.a. diffusion

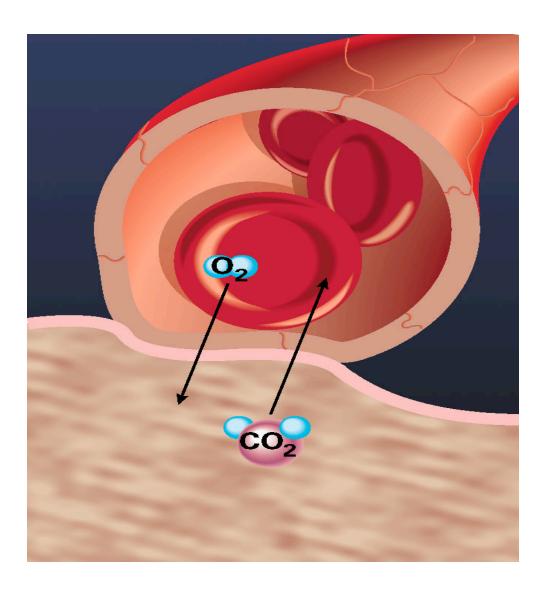


THE PHYSIOLOGY OF PERFUSION

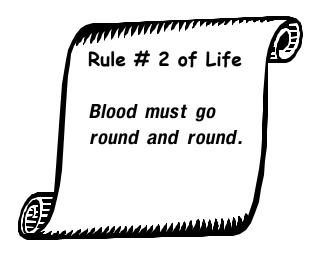

Fick principle: Oxygen Transport
In order from adequate cellular perfusion to occur, the following must be present:

- Adequate number of Red Blood Cells (RBC's)
 - ♦ Hemoglobin on the RBC's carry the oxygen molecules
- Adequate O₂
 - ◆ Patient must have adequate O₂ coming in. See Rule of Life #1
- RBC's must be able to offload and take on O_2
 - Some conditions such as carbon monoxide poisoning and cyanide poisoning affect the RBC's ability to bind and release O_2 molecules.
- Adequate blood pressure to push cells

PHYSIOLOGIC BALANCE

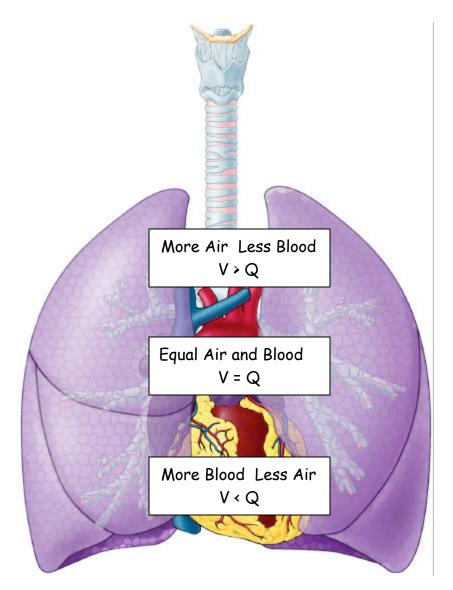

Pathological Conditions

Normal ventilation, poor perfusion: P.E., Arrest

Abnormal ventilation, good perfusion: obstruction, O.P.D., drug OD

Bad ventilation and perfusion: Arrest

Bad exchange area: CHF

Critical Thinking Cases - Designed to illustrate the pathophysiology

Normal ventilation/normal perfusion

Normal ventilation/compromised perfusion

Compromised ventilation/normal perfusion

Compromised ventilation and perfusion

- 1. 26 year-old female patient took an overdose of Valium. She is UNCONSCIOUS. V/S are 110/70, pulse is 64, RR is 12 and very shallow, skin is warm and dry.
- 2. 65 year-old male patient complaining of a sudden onset of right sided chest pain and dyspnea. He has no medical history except for a hip replacement surgery about 3 weeks ago. His lung sounds are clear. VS are B/P 140/78, pulse is 110, and RR is about 20 and normal depth.
- 3. 80 year-old man complains of a sudden onset of severe headache. He has flushed skin, and has obvious facial droop to the left side. He has a history of high blood pressure. V/S are B/P 180/110, pulse is 100 and RR is 16 and normal depth.
- 4. 37 year-old female that was involved in a head on collision. Windshield is starred and the steering wheel is broke. Bruising and crepitus found over the left chest. Pt is unconscious, difficult to bag with absent lungs sounds on the left side. Blood pressure is 60/40; pulse is 130 and weak at the carotid. There is obvious JVD. Skin is cool and clammy.

MODULE TWO: TECHNOLOGY OF CAPNOGRAPHY

The Role of CO2

CO2 is the "Gas of Life"

Produced as a normal by-product of metabolism.

Measurement of EtCO2 (Capnometry)

Qualitative

Color change assay


(CO₂ turns the sensor from purple to yellow)

Quantitative
Gives you a value (EtCO₂)
Respiratory Rate

EtCO₂ **33**mmHg Resp. **13**/min

Waveform Capnography

Features quantitative value and waveform Capnography includes Capnometry

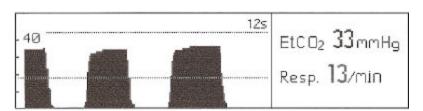

Street Wisdom: "End Tidal CO₂ reading without a waveform is like a heart rate without an ECG recording."

Infrared (IR) Spectroscopy:

Most often used

Infrared light is used to expose the sample

IR sensors detect the absorbed light and calculate a value

Broad spectrum IR beams can be absorbed N_2O and High O_2 levels

Side stream sampling

"First generation devices"

Draws large sample into machine from the line

Can be used on intubated and non-intubated patients with
a nasal cannula attachment

"Second generation devices"
Airway mounted sensors
Generally for intubated patients

Microstream™ Technology

Position independent adaptors Moisture, secretion, and contaminant handling in three ways

- Samples taken from center of line, and in 1/20th the volume
- · Vapor permeable tubing
- Sub micron-multi-surface filters


Expensive parts are protected

Microbeam IR sensor is CO₂ specific

Suitable for adult and pediatric environments.

MODULE THREE: CLINICAL APPLICATIONS OF CAPNOGRAPHY

THE NORMAL CAPNOGRAM

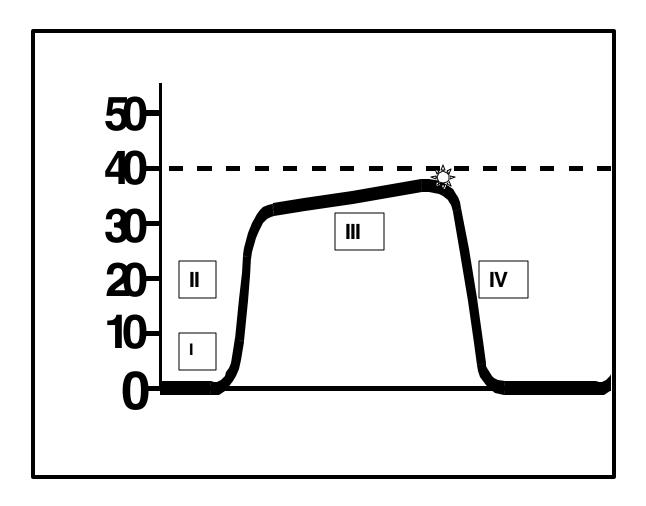

Phase I: Respiratory Baseline

Phase II: Expiratory Upstroke

Phase III: Expiratory Plateau

ETCO2: Peak EtCO2 level

Phase IV: Inspiratory Downstroke

Systematic Approach to Waveform Interpretation

- 1. Is CO₂ present? (waveform present)
- 2. Look at the respiratory baseline. Is there rebreathing?
- 3. Expiratory Upstroke: Steep, sloping, or prolonged?
- 4. Expiratory (alveolar) Plateau: Flat, prolonged, notched, or sloping?
- 5. Inspiratory Downstroke: Sleep, sloping, or prolonged?
- 6. Read the EtCO2
- 7. If ABG is available, compare the EtCO2 with PACO2
 - a. If they are within 5mm/hg of each other then the problem is ventilatory and not perfusion.
 - b. EtCO2 can be used in many cases in lieu of ABG's

The ABC's of Waveform Interpretation!

- A Airway: Look for signs of obstructed airway (steep, upsloping expiratory plateau)
- B Breathing: Look at $EtCO_2$ reading. Look for waveforms, and elevated respiratory baseline.
- C Circulation: Look at trends, long and short term for increases or decreases in Et CO_2 readings


Street wisdom: A patient complains of having difficulty breathing. The pulse oximeter shows 98% on 15 lpm O_2 . As you attempt to listen to lung sounds, they are hard to make out in the back of the ambulance. What benefit, if any could capnography make in the diagnosis and management of this patient?

What is the difference between pulse oximetry and capnography?

SpO₂ = Pulse oximetry - measures oxygenation

EtCO₂ = Capnography - measures ventilation

NORMAL CAPNOGRAPHY

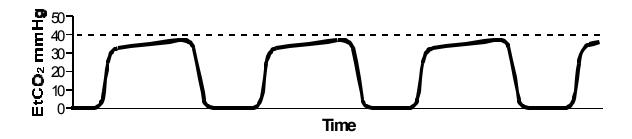


This is a normal capnogram that has all of the phases that are easily appreciated. Note the gradual upslope and alveolar "Plateau"


Point for thought:

List the things a normal capnogram tells you and the things that it does not tell you.

ABNORMAL CAPNOGRAPHY Hyperventilation

This capnogram starts slow and has an $EtCO_2$ reading that is normal. Notice as the rate gets faster, the waveform gets narrower and there is a steady decrease in the $EtCO_2$ to below 30mm/hg. Causes of this type of waveform include:

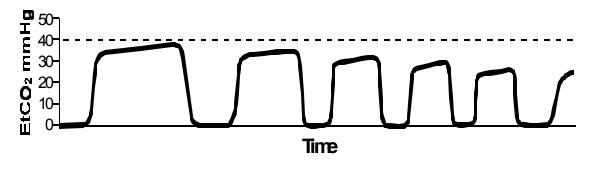

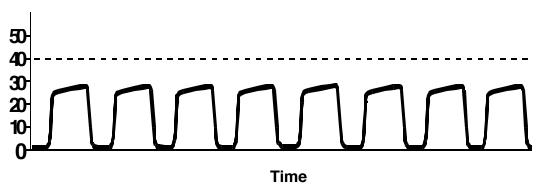


Hyperventilation syndrome

Overzealous bagging

Pulmonary embolism

Hypoventilation

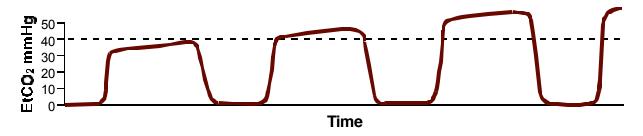


In this capnogram, there is a gradual increase in the $EtCO_2$. Obstruction is not apparent. Causes of this may include:

Respiratory depression for any reason

Narcotic overdose

CNS dysfunction

Heavy sedation

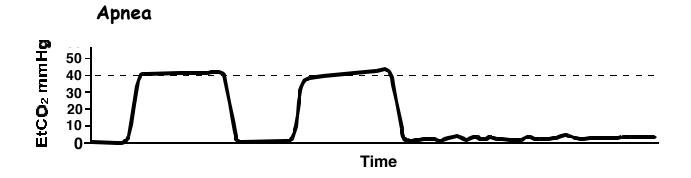

This capnogram shows a complete loss of waveform indicating no CO_2 present. Capnography allows for instantaneous recognition of this potentially fatal condition. Since this occurred suddenly, consider the following causes:


Dislodged ET Tube

Total obstruction of ET Tube

Respiratory arrest in the non-intubated patient

Equipment malfunction (If the patient is still breathing) Check all connections and sampling chambers

Loss of Alveolar Plateau

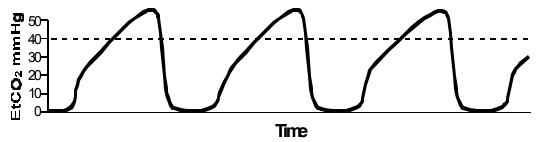

This capnogram displays an abnormal loss of alveolar plateau meaning incomplete or obstructed exhalation. Note the "Shark's fin" pattern. This pattern is found in the following Bronchoconstriction

Asthma

COPD

Incomplete airway obstruction

Upper airway


Tube kinked or obstructed by mucous

Poor perfusion (cardiac arrest)



Time: 30 minute trend

The capnogram can indicate perfusion during CPR and effectiveness of resuscitation efforts. Note the trough in the center of the capnogram. During this time, there was a change in personnel doing CPR. The fatigue of the first rescuer was demonstrated when the second rescuer took over compressions.

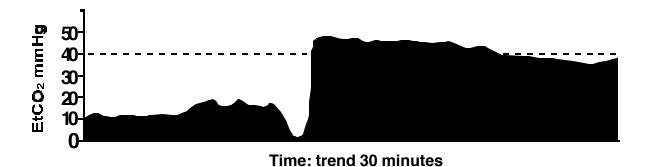


This patient was defibrillated successfully with a return of spontaneous pulse.

Notice the dramatic change in the EtCO2 when pulses were restored.

Studies have shown that consistently low readings (less than 10mm) during resuscitation reflect a poor outcome and futile resuscitation

Elevated Baseline

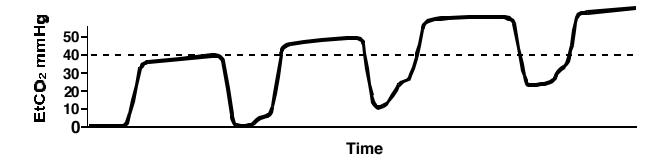

This capnogram demonstrates an elevation to the baseline. This indicates incomplete inhalation and or exhalation. CO_2 does not get completely washed out on inhalation. Possible causes for this include:

Air trapping (as in asthma or COPD)

CO₂ rebreathing (ventilator circuit problem)

What other condition(s) might produce this type of waveform?

Field Clinical Applications for Capnography

Closed Head Injury

Increased intracranial pressure (ICP) tied to increased

blood flow following injury (swelling)

Hypoxic cells produce CO_2 in the brain

 CO_2 causes vasodilation and more blood fills the cranium, increasing pressure.

Hyperventilation is no longer recommended

Ventilation should be geared towards controlling \mathcal{O}_2 levels but not overdoing it.

Obstructive Pulmonary Diseases

Asthma, COPD

Waveform can indicate bronchoconstriction where wheezes might not have been heard

Monitor the effectiveness of bronchodilator therapy

Tube Conformation

Capnography will detect the presence of CO_2 in expired air conforming ETT placement

No longer acceptable to use only lungs sounds to confirm A dislodged tube will be detected immediately with capnography

Kinking or clotting tubes can also be detected In cases of ventilator use, capnography can detect problems in rebreathing.

Perfusion

Capnography can be set up to trend $EtCO_2$ to detect the presence or absence of perfusion

Is a proven predictor of those who do not survive resuscitation

When an ABG is available, can detect ventilatory or perfusion problems

MODULE FOUR: CASE SIMULATIONS AND EVALUATION

CASE # 1


Presentation

Patient is a 65 year old male complaining of crushing substernal chest pain. He rates the pain as a 10 on a scale of one to ten. He denies and shortness of breath or any other complaints. He has a history of cardiac disease and asthma.

Clinical Situation

V/S: 130/80, Pulse is 100, RR is about 20 SpO_2 is 96%, Et CO_2 is 40

Cardiac Monitor shows Sinus Tachycardia His capnogram is as follows.

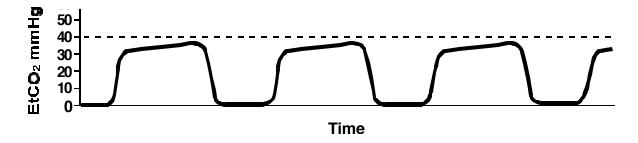

Questions:

Is the EtCO₂ within normal limits?

Is the waveform normal or abnormal? Why or Why Not?

What can you deduce about the ventilation status?

Α

В

 $\boldsymbol{\mathcal{C}}$


CASE #2

Presentation

Patient is a 25-year-old male patient with a history of asthma. He has been compliant with his medications until he ran out of albuterol. Today, while at a basketball game, he suddenly gets short of breath. He does not have his albuterol inhaler with him. He presents sitting in the bleachers, in minor respiratory distress. It is noisy and hard to hear lung sounds.

Clinical Situation

B/P 120/76Pulse 100RR - 14Sp O_2 94 Et CO_2 is 50

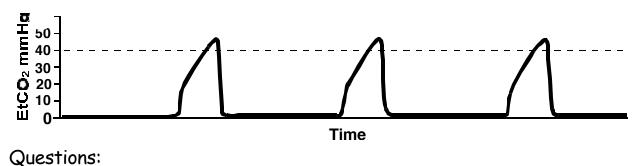

Questions

Is the EtCO₂ within normal limits?

Is the waveform normal or abnormal? Why or Why Not?

What can you deduce about the ventilation status?

Α

В


C

CASE #3

Presentation

You and your partner are working a cardiac arrest and are successful in resuscitation. The Patient is still unstable and the decision is made to load and go because of the very short transport time to the ED. He is intubated and $\rm EtCO_2$ confirmed with good waveform and an $\rm EtCO_2$ of about $\rm 42mm/hg$. The patient is not breathing on their own.

Clinical Situation:

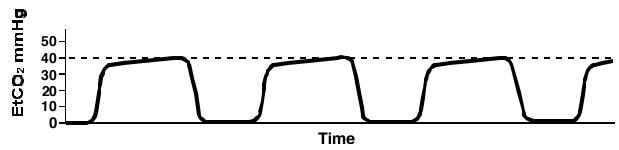

B/P is 100/70


Pulse is 88

RR assisted

SpO₂ is 100% on 15lpm via NRB mask

EtCO2 is 40-42

After loading him into the ambulance, the first responders resume ventilation. The capnography alarm sounds and the following waveform is seen:

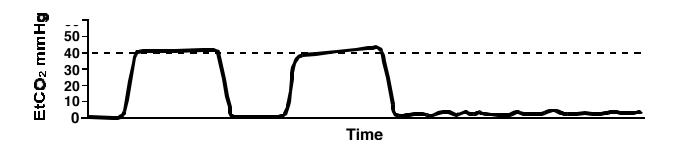

CASE # 3 Continued

Questions:

Is the $EtCO_2$ within normal limits?

Is the waveform normal or abnormal? Why or Why Not?

What can you deduce about the ventilation status?

Α

В

C


CASE # 4

Presentation

You have a 30-year-old female who was in status seizures. Your partner administers Valium to halt the seizures. The patient appears to be post-ictal but is slow to respond fully.

Clinical Situation

B/P is 114/68Pulse is 96 RR is 12 SpO₂ is 98 on 6 lpm nasal cannula Glucose is 100 EtCO₂ is as follows

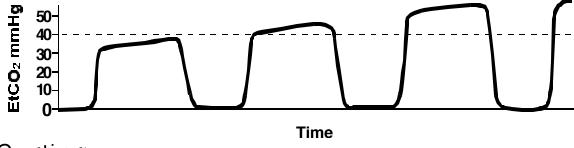

Questions:

Is the EtCO2 within normal limits?

Is the waveform normal or abnormal? Why or Why Not?

What can you deduce about the ventilation status?

A

В

C


CASE # 5

Presentation

It's 3 am and you are called to a residence for a 60 year old man that is in respiratory distress. You find the gentleman sitting up on his bed with feet dangling off the end. He presents in obvious distress and cannot speak words due to the distress. His lung fields are very diminished with crackles heard. He is pale and diaphoretic and appears to be getting weaker. Family members tell you that he has a bad heart and takes a "heart pill", and a "water pill". Pt becomes obtunded with labored breathing. They still have a gag reflex.


Clinical Presentation

BP is 158/90 HR is 130 RR is labored SpO₂ is 88% EtCO₂ as follows

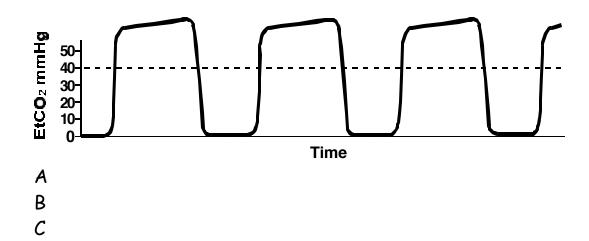

Case #5 Continued

The decision is made to nasally intubate this patient. The tube is passed although the lung sounds are so diminished they are hard to hear. The Pulse Ox offers no change, however, the capnogram shows the following:

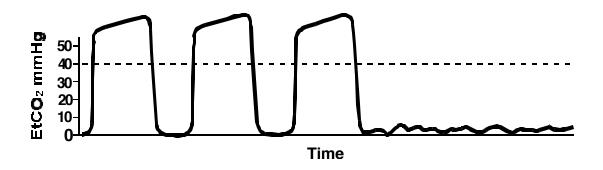

Questions: Scenario 5

Is the $EtCO_2$ within normal limits?

Is the waveform normal or abnormal? Why or Why Not?

What can you deduce about the ventilation status?

A

В

C

Closing Remarks,From one paramedic to another....

Capnography represents another great stride in the advances in technology and medicine that have made way into the field. Not since the cardiac monitor and paramedics manually reading ECG strips has one device had the ability to benefit such a wide variety of patients.

For years, Anesthesiologists have used waveform capnography as their standard for monitoring the vital functions of patients. Now, the technology allows a smaller version to be used by paramedics.

And now, YOU are ready to do this! Think of the incredible difference this can make in the care of your patients.

To summarize, why do you need waveform capnography?

Ventilation Vital Sign

Confirmation of tube placement

Constant monitoring of airway, ventilation and perfusion

Bronchoconstriction in Obstructive airway disease

Any respiratory patient

Closed head injury to guide the careful elimination of CO2

Progressive monitoring of perfusion and ventilation

Why a color change device isn't enough

Only confirms the presence of CO2 not the amount

Can't monitor the patient

Why a quantitative device is not enough

While a number is better than just a color change,

It can't detect bronchoconstriction

It can't trend the level of CO2

There are many different brands and technology out there. I have tried to present a fair and unbiased account of my extensive research into the subject. Take each device for a "test drive" and make your mind up.